

Immunological Test Kit With Immunologically Invisible Carrier

By Bo QIU, Michael BRUNNER, Leonard SIGAL,
Guobao ZHANG, Michael KATZ and Stanley STEIN

Cross References

5 This application claims priority from Stanley STEIN
et al., "Highly Sensitive and Specific IgM-Capture...,"
provisional patent filing serial no. 60/242,819, filed 24 Oct.
2000 and Bo QUI et al., "Multiple Epitopes Connected By A
Carrier," Serial No. 09/_____, filed ____ Oct. 2001. The
10 contents of these, together with Bo QUI, "Studies on Polymers"
(unpublished) and Bo QIU et al., "Selection of Continuous
Epitope Sequences," 55 Biopolymers 319 (2001) are incorporated
here by reference.

Government Rights

15 There are no Federal rights in this invention.

Background

Current technology enables correct diagnosis of
certain infectious diseases only after the disease has
progressed to a certain maturity. By that time, however,
20 treatment is more difficult. We have found a way to make
disease diagnosis, even at an early stage, much more sensitive.

Summary

Our invention entails presenting an immunologically
reactive substance (e.g., epitope polypeptide) in multiple
25 copies conjugated to an immunologically invisible carrier.

This basic conjugate has a variety of versions or
embodiments. For example, while we do not prefer it, the
epitope can be substituted or supplemented with any
immunologically reactive substance such as an epitope, antigen
30 (e.g., a polypeptide or nucleic acid) or antibody. Similarly,
we prefer the carrier also connect a reporter moiety to make
detection of the conjugate simpler.


The conjugate so made may then be used in a variety
of ways. For example, we have shown it effective as part of an
35 immunological assay. Alternatively, the conjugate may be used
as a vaccine. Alternatively, the conjugate may be used as an *in*
vivo therapeutic.

00982265.101701

Thus, our basic idea can be used to make, for example, an immunological test kit. The term "immunological test kit" means a test kit which uses immune (e.g., antibody-epitope or antibody-antigen) interaction to test for the presence or absence of an analyte. Currently-known examples include ELISA, capillary immuno-chromatography and column immuno-chromatography. In making an immunological test kit, it may be desirable to conjugate a reporter moiety on the immunologically invisible carrier (e.g., polyethylene glycol). As another example, our basic idea can be used to conjugate several immunologically reactive substances (either several copies of the same substance, or copies of each of several different substances) together using an immunologically invisible carrier, which conjugate can be then used in an immunological test kit.

The immunologically reactive substance(s) can be one or more of the *Borellia burgdorferii* epitope polypeptides we discovered: VQEGVQQEGAQQP-(beta-A)(beta-,4)C; EIAAKAIGKKIHQNGG-(beta-A)(beta-A)C; ISTLIKQKLDGLKNE-(beta-A)(beta-A)C; PWAESPKKPE-(beta-A)(beta-A)C; DKKAINLDKAQQKLD-(beta-A)(beta-A)C; ITKGKSQKSLGD-(beta-A)(beta-A)C; and GMTFRAQEGAFLTG-(beta-A)(beta-A)C. Alternatively, one could use as antigen the nucleic acid coding for one or more of these epitopes. Using such an epitope enables one to make an apparatus for isolating anti-*Borellia burgdorferi* antibody (i.e., a Lyme disease test kit), a vaccine, or a therapeutic. Similarly, the nucleic acid sequences coding for these polypeptides may be useful as antigen, or to make large quantity of polypeptide.

Our basic idea can be made using, as an immunologically invisible carrier, a polyethylene glycol copolymer that we invented. It has the structure:

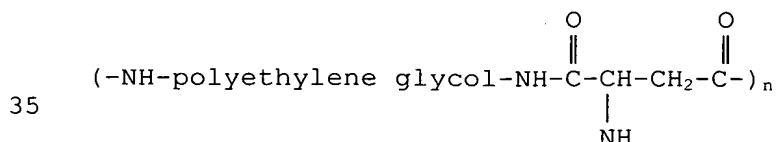

$$\begin{array}{c}
 \text{O} \qquad \qquad \text{O} \\
 \parallel \qquad \qquad \parallel \\
 5 \quad (-\text{NH}-\text{polyethylene glycol}-\text{NH}-\text{C}-\text{CH}-\text{CH}_2-\text{C}-)_n \\
 \qquad \qquad \qquad | \\
 \qquad \qquad \qquad \text{NH} \\
 \qquad \qquad \qquad | \\
 \qquad \qquad \qquad \text{C}=\text{O} \\
 10 \qquad \qquad \qquad | \\
 \qquad \qquad \text{polyethylene glycol} \\
 \qquad \qquad \qquad | \\
 \qquad \qquad \text{O}=\text{S}=\text{O} \\
 \qquad \qquad \qquad | \\
 15 \qquad \qquad \text{CH}_2 \\
 \qquad \qquad \qquad | \\
 \qquad \qquad \text{CH}_2 \\
 \qquad \qquad \qquad | \\
 \qquad \qquad \text{R}
 \end{array}$$

Immunologically Reactive Substance

The success of a specific and sensitive immunoassay largely depends on the strength of antigen-antibody binding and the stability of the complex formed between the antigen and the antibody. The strength of antigen-antibody binding is measured by *affinity*, an intrinsic property of an antigen for a given antibody. To select an epitope peptide is to identify a peptide sequence with high affinity that can bind strongly with specific antibodies.

The stability of complex between antigen and antibody is measured by *avidity*, which is determined by three factors, the *intrinsic affinity* of the antibody for the antigen, the *valence* of the antibody and antigen, and the *geometric arrangement* of the interacting components. Thus, our invention works best when affinity, avidity and specificity (e.g., cross-reactivity) are used to first select an appropriate epitope(s). After the specific epitopes are selected, they can be made as desired (e.g., purified from natural protein or synthesized).

The sensitivity of an immunoassay relies on providing enough of each epitope and on having the right orientation and conformation of the epitope. Thus, we prefer the epitope peptides be modified as necessary to assume the right orientation and conformation to obtain a strong antigen-antibody binding.

Whole antigen or antibody may be used instead of epitope, to mount to the carrier molecule. If mounting antibody on the carrier, the antibody-carrier complex can be used to trap antigen or epitope analyte in the test solution.

Multiple Copies

Epitopes are specific, but have a key shortcoming. The affinity of epitope peptides to anti-protein antibodies can be 100 to 1,000 times weaker than that of the whole antigen (whole protein). Thus, the affinity between a single epitope and the serum antibody might not be strong enough to endure the vigorous washing steps in an immunoassay.

To address this problem, we use multiple copies of each epitope, connected together with a "carrier." Connecting multiple copies of epitope peptides enable the epitopes to form multivalent interactions between two or more Fab fragments of the antibody. This creates a synergistically greater binding strength. More specifically, binding strength increases, perhaps exponentially, with the number of additional copies of epitope connected to the carrier.

For example, an epitope alone may have an antibody affinity 100 times weaker than the native antigen. The same

09982265.101701

epitope, however, if provided in pairs (*i.e.*, two copies of the epitope connected together), might have affinity only 10 times weaker than the native antigen. Further, the same epitope provided in trios (*i.e.*, three copies of the epitope connected together) might have native-strength affinity. We believe this effect especially true where the target antibody is IgM, itself a pentamer.

Immunologically Invisible Carrier

We call the material that connects the various copies of the epitope a "carrier" molecule. Any molecule that can bind more than one copy of an epitope can function as a "carrier." Examples include keyhole limpet hemacyanin, albumins such as serum albumin (*e.g.*, bovine serum albumin, mouse serum albumin, rabbit serum albumin) and ovalbumin, and polyethylene glycol derivatives. These materials can each bind multiple copies of an epitope.

Of these carriers, however, most are unsuitable because they are immunologically "visible," that is to say, they react in an immunological test (even without epitope present) to create a statistically significant increase in (sometimes random) background reactivity. Albumin and limpet hemacyanin tend to stick to ELISA plates. Thus, when using these proteins as carriers, the carrier itself adheres to the ELISA plate in quantity sufficient to cause an elevated background. This problem is particularly significant in developing diagnostic assays for disease where the serum antibody level is relatively low and the signals thus barely detectable. The elevated background compromises the signals, ruining the assay sensitivity and specificity.

Our invention is thus limited to "immunologically invisible" carriers. Excluded from the term "immunologically invisible" are full length albumins and keyhole limpet hemacyanin, because these are not immunologically "invisible."

Biocompatible Polymers

Immunologically invisible carriers are carriers which do not generate statistically significant background

0982265-101701

immunological reactivity. Immunologically invisible carriers include, for example, biocompatible polymers.

Such polymers are known in the art. General reviews of such compounds include Langer, R., "Biomaterials in Drug Delivery," 33 ACC.CHEM.RES. 94 (2000); and Langer, R., "Tissue Engineering," 1 MOL.THER. 12 (2000). One example of such an immunologically invisible compound is a N-vinylpyrrolidone-methyl methacrylate co-polymer, perhaps with added polyamide-6. Buron, F. et al., *Biocompatible Osteoconductive Polymer*, 16 CLIN.MATER. 217 (1994). Another example is poly(DL-lactide-co-glycolide) capsules. Isobe, M. et al., *Bone Morphogenic Protein Encapsulated with a Biodegradable and Biocompatible Polymer*, 32 J.BIOMED.MATER.RES. 433 (1996). Another example is a 70:30 ratio mixture of methylmethacrylate:2-hydroxyethyl methacrylate. Bar, F.W. et al., *New Biocompatible Polymer Surface Coating*, 52 J.BIOMED.MATER.RES. 193 (2000). Another example is 2-methacryloyloxyethyl phosphorylcholine, perhaps with polyurethane. Iwasaki, Y. et al., *Semi-Interpenetrating Polymer Networks...*, 52 J.BIOMED.MATER.RES. 701 (2000). Polyvinyl pyrrolidone may also be used, as may polyethylene glycol and its derivatives. Other biocompatible polymers are known in the art. E.g., Haisch, A. et al., *Tissue Engineering of Human Cartilage Tissue*, 44 HNO 624 (1996); Ershov, I.A. et al., *Polymer Biocompatible X-Ray Contract Hydrogel*, 2 MED.TEKH. 37 (1994); Polous, I.M. et al., *Use of A Biocompatible Antimicrobial Polymer Film*, 134 VESTN.KHIR.IM.II GREK. 55 (1985).

In addition to such synthetic polymers, immunologically invisible biological materials may be used. An example is calcium alginate, such as purified high guluronic acid alginates. Becker, T.A. et al., *Calcium Alginate Gel*, 54 J.BIOMED.MATER.RES. 76 (2001). Genetically engineered protein polymers also may be acceptable. Buchko, C.J. et al., *Surface Characterization of Porous, Biocompatible Protein Polymer Thin Films*, 22 BIOMATERIALS 1289 (2001); cf. Raudino, A. et al., *Binding of Lipid Vesicles...*, 231 J.COLLOID.INTERFACE SCI. 66 (2000).

09082265 101701

Such compounds may lack functional groups useful for attaching the desired immunologically reactive substance to the carrier. Thus, it may be desirable to use not the pure polymer, but a co-polymer having appended functional groups. The functional groups may then be filled with the desired immunologically reactive substance.

As immunologically invisible carrier, we prefer polyethylene glycol and its derivatives. We thus now discuss it in some detail.

10 *Polyethylene Glycol*

Polyethylene glycol (often simply called "PEG") is a water soluble, non-immunogenic, biocompatible material. When used as a carrier, the useful properties of polyethylene glycol with respect to the appended moiety include improved solubility, increased circulation lifetime in bloodstream, resistance to proteases and nucleases, etc. The large molecular weight of polyethylene glycol makes it very easy to separate the final conjugates from excess epitope peptide and other small-size impurities. Polyethylene glycol does not aggregate, degrade or denature. Polyethylene glycol conjugates are thus stable and convenient for use in diagnostic assays.

While the polyether backbone of polyethylene glycol is chemically inert, the primary hydroxyl groups on both ends are reactive and can be utilized directly to attach immunologically reactive substances. These hydroxyl groups have been transformed into more reactive functional groups for conjugation purposes. Such polyethylene glycol derivatives possess only two functional groups on the ends. This limits the number of conjugations to just two. We thus prefer a polyethylene glycol derived polymer system with multiple functional groups for epitope peptide attachment.

We made a new polyethylene glycol with multiple functional groups and a favorable geometric arrangement to achieve strong and stable antigen-antibody binding for the selected epitope peptides. We used α,ω -diamino-polyethylene glycol to copolymerize with amino group-protected aspartic acid

09982265.101701

to obtain a new polyethylene glycol-aspartic acid copolymer. Multiple attachment sites become available for conjugation through the pendant amino groups of the aspartic acid residue upon removal of the protection (Figure 1).

5 To allow the attached epitope peptides to assume a favorable geometric arrangement for antibody binding, we used a long arm cross-linker for attaching the epitope peptides to the amino groups, so that the attached epitope peptides can be positioned far enough from the polymer backbone to avoid steric
10 hindrance. We used a heterobifunctional polyethylene glycol-based cross-linker, NHS-polyethylene glycol-VS, as the cross-linker for epitope peptide conjugation.

The conjugation of epitope peptides may use thiol-specific chemistry under mild conditions. The easiest strategy
15 for peptide conjugation is to add an extra amino acid on either the amino or carboxyl terminus of the peptide to allow one-site coupling to the carrier. In our study design, a cysteine residue, followed by two β -alanine residues, was incorporated at the C-terminus of each epitope peptide during solid phase
20 peptide synthesis. Putting two more β -alanine residues between the conjugation anchor, cysteine, and the epitope peptide is used as a precaution to generate further flexibility of the linear peptides, and therefore help them to adopt the optimal conformations for stronger antibody binding. The N-terminus of
25 the peptides needs to be capped in order to remove charges associated with free amino groups and thereby mimicking the real environment in the protein.

To conjugate epitope peptides to the polymer backbone, a two step approach can be used. A heterobifunctional
30 cross-linker, NHS-polyethylene glycol-VS can first react with the reserved amino groups in the reporter-labeled polymer carrier through, the NHS groups. After purification to remove excess cross-linker, cysteine-containing epitope peptides can then react readily with vinylsulfone groups (VS) to complete the
35 conjugation. The final polyethylene glycol -peptide conjugates containing multiple copies of epitope peptides and several

09982265.101701

Reporter

A commonly used reporter molecule in immunoassay is biotin. Its corresponding N-hydroxysuccinimide ester (NHS) with extended spacer is chosen for our carrier-peptide conjugate preparation. We did this because the NHS group can react readily with the pendant amino groups of the polyethylene glycol-aspartic acid copolymer under mild conditions. The extended spacer arm can help lower steric hindrance and thus facilitate assay detection. Since biotin detection system is extremely sensitive, a few label molecules should suffice to give satisfactory signals. Therefore, only a small portion of attachment sites in the carrier is needed to attach reporter molecules so that a large portion of the attachment sites can be reserved for the epitope peptides to generate polyvalent antigen with improved antibody binding and to improve the sensitivity of the immunoassay.

Alternatively, the reporter molecule can be put on the N-terminus of the epitope peptides during the solid phase peptide synthesis. The reporter molecules can thus serve as the capping groups of the peptides and as the reporter groups of the conjugates simultaneously. By putting the reporter groups both on the polymer backbone and on the epitope peptides, the assay signal can be further enhanced (Figure 3). Care must be taken to not block the epitope from contacting and binding to the antibody. Multiple copies of the reporter groups attached to the carrier amplify the assay signal. Other reporters or labels

(e.g., colloidal metal, carbon black, latex beads) are known in the art and may alternatively be used.

Uses

Once made, our carrier-epitope conjugates can be used for a variety of things. For example, our conjugates can be used in immuno-chromatography, the specific kind of chromatography selected depending on one's goals. Column chromatography, for example, can be done with our conjugates used to isolate and purify a desired antibody in quantity. Alternatively, capillary chromatography can be done with our conjugates, to detect low levels of antibody in a sample. Similarly, ELISA can be done with our conjugates, to detect low levels of antibody in a clinical sample. We actually used our conjugates to make such an immunodiagnostic kit, so we will now discuss how to make such a kit in some detail.

Detailed Description of our Preferred Embodiment

Our preferred embodiment of our invention entails four parts: 1) the selection of specific epitopes by epitope mapping; 2) the design and synthesis of a carrier molecule with multiple attachment sites; 3) the preparation of multivalent carrier-peptide conjugates with one or more reporter groups; and 4) the use of the prepared carrier-peptide-reporter conjugates in an immunological assay. Here is how you can use of our preferred embodiment to make an indirect IgM-capture ELISA effective for the diagnosis of Lyme disease at its earliest stage.

Epitope Mapping by SPOTS

All concurrent peptide sequences were generated using computer software provided by the manufacturer (Genesys) with the SPOTS kit. By providing a protein sequence, desired length of each peptide and offset of amino acids for each peptide, the program can edit peptide sequences to be assembled on SPOTS membrane and provide the amino acid addition schedule for each synthesis cycle.

0992265-101701

To start the peptide synthesis on the membrane, pre-weighed Fmoc-amino acid active esters were dissolved in DMF and pipetted to appropriate spots on the membrane based on the generated synthesis schedule. Double coupling was done for each cycle to ensure the completion of the reaction. All the Fmoc-amino acid active esters, except Arginine, are relatively stable and can be dissolved in DMF for use of several cycles in the same working day, as long as they are stored at -20 °C between each addition. Due to its intrinsic instability, the Fmoc-Arginine active ester must be dissolved just before use and a fresh aliquot must be used for each coupling cycle. The initial color of all spots on the membrane was blue which is produced by bromophenol blue in the presence of the free amino groups on the de-protected amino acids.

As coupling proceeds with the addition of Fmoc-amino acid active esters, the spots change to different colors for different amino acids. For example, Asparagine and Threonine change to green, Serine changes to yellow. The color change can be regarded as a sign that the coupling is taking place. After coupling an amino acid the membrane was washed 3x20 mL DMF for 2 minutes each time to remove excess active esters.

Then, acetic anhydride was added to acetylate any uncoupled amino groups to ensure no formation of deletion sequences. As all free amino groups are capped by acetylation, the remaining blue color disappeared. The membrane was washed 3x20 mL DMF and then 20 mL of 20% piperidine in DMF was added to remove Fmoc protecting groups. After washing membrane 5x20 mL DMF, 200 µL of 1% bromophenol blue solution was added to 20 mL DMF and this solution was added on the membrane. Due to piperidine removal of the Fmoc groups, the spots turned blue leaving the surrounding membrane white and the solution yellow. The membrane was washed 3x20 mL with methanol. After air drying on a sheet of folded filter paper, the membrane is ready for the next coupling cycle. This procedure was repeated for all but the final coupling cycle of the synthesis.

For the final cycle, piperidine treatment was carried

09982265.101701

09982265 101701

out right after the double coupling of active esters and DMF washing. Bromophenol blue solution was then added to obtain blue color for all spots and finally the peptides on each spot were capped by acetylation. After synthesis and acetylation, the protecting groups present on the side chains of the amino acids must be removed. For side chain deprotection, 5 mL of DCM was mixed with 5 mL TFA. The mixed solution was added immediately onto the air-dried membrane and the cleavage reaction was allowed to proceed for 1 hour. The membrane was then washed with 3x20 mL DCM, 3x20 mL DMF, and 3x20 mL methanol. The membrane was air-dried and stored in a sealed plastic bag in the freezer (-20 °C) until required for SPOTS analysis.

For analysis, the SPOTS membrane was first blocked with 20 mL of TBS-blocking buffer overnight at room temperature. The membrane was washed with 20 mL; Tris buffered saline (TBS) containing 0.05% Tween-20 (T-TBS). The serum sample (Lyme disease or control) was diluted in 20 mL TBS-blocking buffer to 1:100. This diluted test antibody solution was added to the membrane and rocked for 3-4 hours at room temperature. The membrane was washed with 3x20 mL T-TBS for 10 minutes each wash. Then, 100 µL of P-galactosidase conjugated anti-human (G+M+A) secondary antibody was diluted with 20 mL of TBS-blocking buffer. This was added to the membrane and rocked for 2 hours at room temperature.

During this time, the signal development solution was prepared as follows: Dissolve 4.9 mg BCIG in 100 µL DMF and 100 mg potassium ferricyanide in 1 mL MilliQ water. Add BCIG solution and 100 µL of potassium fericyanide solution into 10 mL of phosphate buffered saline (PBS) containing 10 µL of 1 M magnesium chloride solution. After the incubation of the secondary antibody solution, wash the membrane 2x20 mL T-TBS followed by 2x20 mL PBS, then add the prepared signal development solution to the membrane and rock at room temperature until blue spots appear. Allow the color to develop for 40 to 50 minutes until a point at which there is a clear distinction between positive and negative spots. Pour off the

signal development solution and wash the membrane with 2x20 mL PBS. Photograph the stained membrane to provide a permanent record.

The SPOTS membrane must be regenerated after analysis of each serum sample to remove I bound proteins before storage or re-probing. To regenerate the membrane, it was washed with 5x20 mL MilliQ water and then 3x20 mL DMF followed by another 2x20 mL MilliQ water. Then, 20 mL, of regeneration buffer A (485.0 g urea, 10.0 g SDS and 1 mL 2-mercaptoetbanol in 1 L of MilliQ water) was added and the membrane was incubated for 10 minutes at room temperature. The process was repeated twice. Then 20 mL of regeneration buffer B (Mix 400 mL of MilliQ water and 500 mL ethanol, add 100 mL of acetic acid to above solution) was added and the membrane was incubated for 10 minutes at room temperature. The process was repeated twice. Finally, the membrane was washed with 2x20 mL methanol and air-dried. The membrane was stored in a sealed plastic bag in the freezer (-20 °C) until the next analysis.

Synthesis, Purification and
Characterization of Epitope
Peptides

All 7 epitope peptides (Table 1) were synthesized manually on PAL™ resin (0.34 mmol/g, 0.1-0.2 mmol scale) in a polypropylene column (Bio-Rad Laboratories, Hercules, CA). DMF (3 ml) was added to swell the resin for 20 min. After Fmoc de-protection with 20% piperidine in DMF for 2x20min, the resin was rinsed with 3 ml of DMF three times, 3 ml of methanol three times, then dried in air. The coupling was achieved by adding 3-fold molar excess of each amino acid, mixed with equimolar amounts of BOP and HOBt in 3 ml of DMF containing 1% (v/v) DHEA. Coupling proceeded at room temperature for 4 hours.

After coupling, the resin was washed with DMF and methanol and air-dried. A sample of the resin was tested with Kaiser ninhydrin reagent (1: 1: 1 v/v/v 0.2 mM KCN in pyridine, 4 mg/ml of phenol and 5% ninhydrin in butanol) at 10 °C for 3 min (Kaiser et al., 1970; Sarin et al., 1981). If the resin

09982265-101701

showed blue color, double coupling would be conducted for another 4 hours to drive the reaction to completion. The resin was capped using 4 mL of DMF, 400 μ L of acetic anhydride and 80 μ L of triethylamine for 4 hours to eliminate any un-reacted amino groups.

The coupling procedure was repeated until the desired peptide sequence was obtained. When the assembly of the peptide sequence was complete, the N-terminus of all epitope peptides was capped with long chain biotin to serve two purposes simultaneously. The first purpose is to remove the charge associated with the free amino group of the N-terminus, thus to mimic the real environment in the natural protein sequence. The second purpose is to use the biotin as the detection label for biotin-avidin binding in ELISA.

Table 1
 Synthesized Epitopes

Peptide	Sequence	
FLA, AA 211-223	VQEGVQQEGAQQP-(beta-A)(beta-,4)C	1639.8
OspC2, AA71-86	EIAAKAIGKKIHQNNG-(beta-A)(beta-A)C	2274.3
OspC3, AA 104-118	ISTLIKQKLDGLKNE-(beta-A)(beta-A)C	2282.3
OspC10, AA 198-207	PWAESPCKPE-(beta-A)(beta-A)C	1762.7
P83-1, AA296-310	DKKAINLDKAQQKLD-(beta-A)(beta-A)C	2310.3
P83-3, AA431-442	ITKGKSQKSLGD-(beta-A)(beta-A)C	1843.8
P39, AA129-142	GMTFRAQEGAFITG-(beta-A)(beta-A)C	2067.9

Long chain biotin was selected to reduce any possible conformational hindrance for high-avidity biotin-avidin binding. All peptides were cleaved from the resin with trifluoroacetic acid (TFA) / thioanisole / ethanedithiol (EDT)/anisole (90/5/3/2%, v/v) at 1 mL/100 mg resin for 2 hours at room temperature. The cleavage mixture was filtered through glass wool, which was then rinsed with TFA twice. The filtrates were combined and evaporated under an Argon stream to reduce the volume to about 1-2 mL, then precipitated by adding drop-wise

00982265.101701

into 10 times volume of ice-cooled diethyl ether. The white precipitate was washed with cold diethyl ether five times to remove scavengers. Crude peptides were purified by reverse phase HPLC under acidic condition (0.1% TFA), because cysteine was incorporated in all epitope peptides for conjugation purpose and the availability of free thiol groups in cysteine is critical for conjugating epitope peptides onto PLC copolymer backbone. The acid condition can help to prevent or minimize the oxidation of the free thiol groups. After HPLC purification, the tubes containing the epitope peptides were flushed with Argon stream, capped, wrapped with paraffin, and stored dry in the refrigerator (4 °C). The purified epitope peptides were characterized by amino acid analysis and mass spectrometry.

15 Synthesis and Purification of
polyethylene glycol-Aspartic Acid
Copolymers

Amino group protected L-Aspartic acid (Boc-Asp-OH) (BACHEM, King of Prussia, PA) and α,ω -diamino-PEG (NH₂-PEG-NH₂, Shearwater Polymers, Huntsville, AL) were copolymerized based on carbodiimide reaction in the presence of 4-(dimethyl amino)-pyridine (DMAP) and p-toluenesulfonic acid monohydrate (PTSA) as catalysts. In a typical preparation, NH₂-PEG-NH₂ (680 mg, 2×10^{-4} mol) and Boc-Asp-OH (46.6 mg, 2×10^{-4} mol) were dissolved in 20 mL methylene chloride with stirring. DMAP (12.2 mg, 1×10^{-4} mol) and PTSA (19.0 mg, 1×10^{-4} mol) were added. To this solution 1,3-diisopropylcarbodiimide (DIPC) (15.6 mL, 1×10^{-3} mol) was added at 0 °C under stirring. The reaction flask was sealed with a rubber stopper assembled with an Argon balloon. The reaction was allowed to continue at room temperature with stirring until the reaction mixture became viscous.

The reaction mixture was precipitated in 10 volumes of ice-cold ethyl ether to obtain the white polymer product. The polymer was washed three times with ice-cold ethyl ether and the polymer product was collected by filtration or centrifugation. The polymer was dried under an Argon flow, re-

09982265-101701

dissolved in MilliQ water and purified by dialysis using Spectra/For™ Spectrum cellulose ester membrane (MW 12-14,000 Da) for 24 h. After lyophilization, the polymer was treated with TFA for 3 hours to remove all the Boc protecting groups. The de-protected polymer solution was then precipitated in 10 volumes of ice-cold ethyl ether, washed three times with ice-cold ethyl ether and dried under vacuum. The molecular weight of the resulting PEG copolymer was measured by size exclusion chromatography.

10 Preparation of polyethylene glycol-
Peptide Conjugates

To a solution of PEG copolymer in 50 mM carbonate-bicarbonate buffer (pH = 8.5) was added 0.5 equivalent (relative to the amino groups in the polymer) of NHS-LC-Biotin in DMSO. The mixture was stirred at room temperature under Argon overnight. After about 10 hours of reaction, approximately 30% of the amino groups in the PEG, copolymer were reacted and linked to biotin molecules. A fluorometric assay, using a fluorogenic reagent, Fluram, was employed to check the extent of the biotinylation reaction. In brief, 100 µL of PEG copolymer solution was saved before adding the biotinylation reagent and diluted 10x in 0.2 M borate buffer, pH 8.5) as reference. When reaction was complete, 100 mL of reaction mixture was taken and diluted 10x in 0.2 M borate buffer (pH 8.5) as sample.

For fluorometric assay, 50 mL of Fluram solution (15 mg Fluram dissolved in 25 mL acetonitrile) was added to 150 µL of diluted reference, 150 µL of diluted sample and 150 µL of blank (0.2 M borate buffer, pH 8.5), respectively, in separate wells of a microtiter plate. After mixing immediately by pipetting up and down several times, fluorescence was read on a Fluorescence Multi-Well Plate Reader (CytoFluor™ 11, PerSeptive Biosystems) with the excitation wavelength set at 400 nm and the emission wavelength set at 460 nm. The biotin labeled PEG copolymer was purified by a Pharmacia Superdex-75 column and then reacted with 3 molar equivalents of hetero-bifunctional NHS-PEG-VS (MW 2000 Da), relative to free amino groups remaining

09982265-101701

in biotin-labeled PEG copolymer.

The latter reaction, which was also monitored by the fluorometric assay, was complete after 4 hrs at room temperature (25 °C). The fluorometric assay procedure was similar to that described above. The final fluorescence reading was equal or close to the blank reading, suggesting that (all amino groups in the PEG copolymer had been successfully derivatized. The reaction product was purified through a Pharmacia Superdex-75 column or by membrane dialysis. For peptide conjugation, 5 molar equivalents of peptide relative to the available vinylsulfone (VS) groups in the PEG copolymer were added to the activated polymer solution, and these were allowed to react at 4°C overnight. The final Biotin-PEG-peptide conjugate was purified by the Pharmacia Superdex-75 column or by membrane dialysis, and concentrated to about 1 mg/mL using a Centricon™ ultrafilter (mw 10,000 Da). Aliquots were stored as the stock antigen solution in the freezer (-20 °C) until needed.

The Enzyme-Linked Immuno-Sorbent Assay

ELISA is a simple but very sensitive immunoassay. It involves the following basic steps: An antigen is bound to a solid phase material, usually a 96-well plastic plate. The solution containing the antibody to be detected (usually serum) is added to the well having the immobilized antigens. Unrelated, unbound antibody is then washed away. A second antibody, which is an anti-immunoglobulin antibody linked with an enzyme, is then added to the wells. Then the substrate for the enzyme is added to the above reaction mixture and the amount of enzymatically altered substrate is measured. The enzyme and substrate are chosen so that enzymatic modification of the substrate produces a change in color of the substrate solution. The amount of changed substrate (which may be measured with a spectrophotometer) is proportional to the amount of antibody bound to the immobilized antigen.

There are generally two types of ELISA formats: direct and indirect. In a direct ELISA, antigens first bind to

the well surface of the plates, and then the bound antigens interact with the test antibodies and give the signals. In an indirect ELISA, the plates are first coated with antibodies that can capture antigens. The captured antigens can then interact with the test antibodies and give the signals.

Many modifications of the above basic technique can be used depending on the nature of the sample, availability of reagents and the precision and sensitivity required. For example, one may use a biotinylated antibody followed by enzyme-conjugated avidin or streptavidin. The avidin-biotin method results in an amplified effect since many biotin molecules may be attached to a single second antibody molecule and multiple avidin molecules can then bind subsequently to the second antibody. For this reason, the avidin-biotin method is particularly sensitive.

i) IgM Capture ELISA

In an IgM-capture format, IgM antibodies are captured or bound to the test support, such as an ELISA plate. A representative portion of all IgM antibodies, including disease specific and unrelated IgM antibodies, are captured. All other classes of antibodies are removed.

In a direct-capture test, the antigens are immobilized on the surface of the plate. In an indirect-capture test, the antigens are present in the test solution and interact with the antibodies captured or bound to the ELISA plate.

When the captured IgM antibodies are exposed to the prepared PEG-peptide conjugates, these Lyme disease specific epitope conjugates will only bind to Lyme disease specific IgM antibodies. If no Lyme disease specific IgM antibodies are present, all conjugates will be washed away and no signal can be detected. As a result, a negative result is obtained. Clearly, this indirect IgM capture ELISA format, combined with using the Lyme disease specific conjugates as antigens, increases the sensitivity and the specificity of detecting Lyme disease specific IgM antibodies, on which a highly sensitive and specific immunoassay can be developed (Figure 4).

0982265.101701

09982265-101701

ELISA plates were coated with 100 μ L/well of affinity-purified goat anti-human IgM antibody (10 μ g/mL) in 0.04 M carbonate-bicarbonate buffer, pH 9.6. Plates were slowly rotated on a Titer Plate Shaker (Lab-Line, Melrose Park, IL) for 2 h at room temperature, and kept at 4°C overnight. The plates were washed three times in a plate washer (ELP 3.5, Biotek, Winooski, Vt.) with PBS-B (10 mM phosphate buffered saline, 0.15 M sodium chloride, containing 0.1% BSA), blocked with 300 μ L/well of PBS-B milk (PBS-B containing 5% nonfat dry milk) for 2 h at 37°C. Serum samples were diluted 1:100 in PBS-B milk, added at 100 μ L/well and rotated at 300 rpm for 1 h. The plates were washed four times with PBS-B and incubated for 1 h with 100 μ L/well of Biotin-PEG-peptide conjugates (diluted to various concentrations in PBS-B milk).

During this time, the avidin-biotinylated peroxidase complex (ABC) was formed by adding one drop (50 μ L) of reagent A (avidin DH) and one drop (50 μ L) of reagent B (biotinylated peroxidase) to 5 mL of PBS-BT (PBS-B containing 0.5 M sodium chloride and 0.1% Tween 20). The ABC reagent was vortexed and kept at room temperature for at least 30 minutes before use. After washing the plates four times with PBS-B, 7 mL of PBS-BT was added to the ABC reagent and 100 μ L of the diluted ABC reagent was added to each well. The plate was rotated at 300 rpm for 30 minutes and washed four times with PBS-B on the Biotek plate washer followed by two more manual washes with plain PBS. During the last wash, the two component 3,3',5,5'-tetramethylbenzidine substrate solution (TMB) was prepared at room temperature. Substrate was added at 100 μ L/well with a repeater pipette (Eppendorf Plus/8), the plate was rotated for 10 minutes to develop the color, and the reaction was stopped by adding 100 μ L/well of 1 M phosphoric acid. The plate was then rotated for 2 more minutes to homogenize the color and then read on an ELISA plate reader (Biotek) set for dual wavelengths (450 and 630 nm).

All seven Biotin-PEG-peptide conjugates were tested as antigens in IgM-capture ELISA individually and as in

combination with a panel of samples containing sera from both
Lyme disease patients and healthy subjects. A group of 12
negative control sera. were tested under the same assay
conditions and the average absorbance plus three standard
5 deviations of these control serum samples was used as the
cutoff.

The index number of each serum sample was calculated
as: Index = Absorbance of individual serum/Cutoff. An index
number of 1.0 or above is taken as a. positive and an index
10 number of 0.8 or below is taken as a negative. Any index number
between 0.8 to 1 .0 is taken as equivocal.

ii) Clinical diagnosis
by IgM Capture ELISA

A panel of sera is tested by IgM capture ELISA using
15 either protein-based antigen (*Borrelia burgdorferi* sonicate) or
our peptide-based antigens. The clinical diagnosis results are
listed in Table 2.

The peptide-based ELISA using the combination of
seven PEG-peptide conjugates identified 31 positive samples from
20 33 culture-proven positive samples, resulting in a diagnostic
sensitivity of 94% (percentage of disease samples correctly
diagnosed). The protein-based ELISA using sonicated *Borrelia*
burgdorferi spirochete picked up 23 samples out of 31 tested
positive sera, yielding a diagnostic sensitivity of 74%.
25 Furthermore, the peptide-based ELISA did not yield any false
positive results with the non-Lyme disease samples giving an
essentially 100% of diagnostic specificity, whereas the protein-
based ELISA gave 6 false positives out of 23 negative samples,
or a diagnostic specificity of 74% (percentage of non-disease
30 samples correctly diagnosed). Thus, the peptide-based ELISA
achieved higher sensitivity and specificity than the protein-
based ELISA.

As our design rationale predicted, the defined
epitope peptides should have less tendency than whole proteins
35 to cross-react with sera from patients with other diseases, such
as syphilis. In order to examine this hypothesis further, a

09982265.101701

panel of serum samples from patients with syphilis infection was tested using the combination of PEG-peptide conjugates (Table 3). Indeed, while 13 out of 25 syphilis samples gave cross-reactive results in the protein-ELISA, none of these tested
5 syphilis samples showed cross-reactivity in our peptide-ELISA when corrected by subtracting serum background (no antigen used in ELISA), indicating that all seven epitope peptides defined in this study are Lyme disease specific and do not cross-react with antibodies against the syphilis spirochete.

10

Summary

In our claims, we use the singular to include the plural (i.e., "a" or "an" means "one or more").

The present invention is not to be limited in scope by the specific embodiments disclosed in the examples which are
15 intended as illustrations of a few aspects of the invention and any embodiments which are functionally equivalent are within the scope of this invention. Indeed, various modifications of the invention in addition to those shown and described herein will become apparent to those skilled in the art and are intended to
20 fall within the scope of the invention. Thus, for example, serum antibodies specific for any disease can be analyzed in order to select disease-specific epitope sequences. Peptides corresponding to these epitope sequences are then synthesized and conjugated in several copies to a multivalent PEG carrier
25 molecule, along with a reporter group such as biotin. We thus intend the legal coverage of our patent to be defined not by the scientific examples we include here, but by the legal claims appended here.

30

09982265 "101701

Table 2. Comparison of Lyme disease diagnosis for a panel of serum samples.

No.	Protein ELISA	Peptide ELISA	Clinical diagnosis	No.	Protein ELISA	Peptide ELISA	Clinical diagnosis
MC-2	P	P	P	NC-1	P	N	N
MC-3(6/18)	N	P	P	NC-2	N	N	N
MC-3(6/26)	P	P	P	NC-3	N	N	N
MC-4	P	P	P	NC-4	N	N	N
MC-7	P	P	P	NC-5	N	N	N
MC-8	N	P	P	NC-8	N	N	N
MC-9	P	E	P	NC-9	N	N	N
MC-10	P	P	P	NC-10	N	N	N
MC-14	P	P	P	NC-11	N	N	N
MC-17	P	P	P	NC-14	N	N	N
MC-23	P	P	P	NC-15	N	N	N
MC-33	P	P	P	NC-16	N	N	N
MC-41	P	P	P	NC-A	N	N	N
MC-59	P	P	P	NC-B	N	N	N
MC-62	P	P	P	NC-C	P	N	N
MC-68	N	E	P	NC-D	N	N	N
MC-70	P	P	P	NC-E	P	N	N
MC-71	E	P	P	NC-F	N	N	N
MC-72	P	P	P	NC-G	N	N	N
MC-73	P	P	P	NC-H	N	N	N
MC-74	N	P	P	NC-LT	P	N	N
MC-91	N	P	P	NC-LN	P	N	N
MC-92	P	P	P	NC-LP	P	N	N
MC-93	N	P	P	MC-SC	P	P	P
MC-100	N	P	P	MC-EL	P	P	P
MC-101	P	P	P	MC-AN	P	P	P
MC-JS	P	P	P	MC-MT	ND	P	P
MC-GR	P	P	P	MC-HA	ND	P	P

P, positive; N, negative; E, equivocal, ND, not determined.

09982265-101701

Table 3. Comparison of ELISA results for serum samples from patients with syphilis.

No.	Protein ELISA	Peptide ELISA	Clinical diagnosis
S-1	CR	NCR	Syphilis
S-2	CR	NCR	Syphilis
S-3	CR	NCR	Syphilis
S-4	CR	NCR	Syphilis
S-5	CR	NCR	Syphilis
S-6	CR	NCR	Syphilis
S-7	CR	NCR	Syphilis
S-8	NCR	NCR	Syphilis
S-9	CR	NCR	Syphilis
S-10	NCR	NCR	Syphilis
S-11	NCR	NCR	Syphilis
S-12	NCR	NCR	Syphilis
S-13	CR	NCR	Syphilis
S-14	NCR	NCR	Syphilis
S-15	NCR	NCR	Syphilis
S-16	NCR	NCR	Syphilis
S-17	NCR	NCR	Syphilis
S-18	NCR	NCR	Syphilis
S-19	CR	NCR	Syphilis
S-20	CR	NCR	Syphilis
S-21	CR	NCR	Syphilis
S-22	CR	NCR	Syphilis
S-23	NCR	NCR	Syphilis
S-24	NCR	NCR	Syphilis
S-25	NCR	NCR	Syphilis

CR, cross-reactive; NCR, non-cross-reactive.